

Ordering Information

Drug Info			For Mutil-Drug Test Panel and Cup	For Single Drug Test Cassette
				
Abbr.	Tests	Specimen	Cut-off (ng/ml)	Cut-off (ng/ml)
AMP	Amphetamine	Urine	500/1000	500
BAR	Barbiturates	Urine	300	300
BUP	Buprenorphine	Urine	10	10
BZO	Benzodiazepines	Urine	300	300
CAF	Caffeine	Urine	6000	\
COC	Cocaine	Urine	150/300	300
COD	Codeine	Urine	300	\
COT	Cotinine	Urine	100	\
EDDP	Methadone metabolites (EDDP)	Urine	300	\
ETG	Ethyl Glucuronide	Urine	500/1000	\
FTY	Fentanyl	Urine	20/100	\
K2	Synthetic Marijuana	Urine	25	\
KET	Ketamine	Urine	1000	1000
6-MAM	6- Monoacetylmorphine	Urine	10	\
MCAT	Methcathinone	Urine	300	\
MDMA	Ecstasy	Urine	500	500
MET	Methamphetamine	Urine	500/1000	500
MOP	Morphine	Urine	300	300
MQL	Methaqualone	Urine	300	\
MTD	Methadone	Urine	300	300
OPI	Opiate	Urine	2000	\
OXY	Oxycodone	Urine	100	\
PCP	Phencyclidine	Urine	25	25
PPX	Propoxyphene	Urine	300	\
SOMA	Carisoprodol	Urine	500/1000	\
TCA	Tricyclic Antidepressants	Urine	1000	1000
THC	Marijuana	Urine	50/20	50
TRA	Tramadol	Urine	100	100



Drug of Abuse Screening Test
Rapid, Reliable & Configurable Drug Testing Solutions



Broad Coverage

Covers classic drugs and new psychoactive substances (NPS).

Integrated Design

Combines sample testing and result reading into one streamlined process.

Smart Packaging

Cutting shipping volume by 30-50%, our eco-friendly packaging slashes your shipping costs and carbon footprint.

Accurate & Rapid

CE-marked and ISO 13485 certified; Results in 5 minutes with high sensitivity and specificity.


Cost-Effective Screening

Affordable solution for large-scale or routine drug testing.

Key Facts

- In 2023, 316 million people worldwide used psychoactive drugs, up from 292 million in 2022.
- Synthetic drug markets—including fentanyl, amphetamines, and cocaine—reached record highs, with 25 million cocaine users reported globally.
- Rapid point-of-care urine tests support early detection, intervention, and workplace safety.
- offering timely, accurate drug screening at scale.

(Source: World Drug Report 2024–2025, United Nations Office on Drugs and Crime (UNODC); WHO Global Health Observatory.)

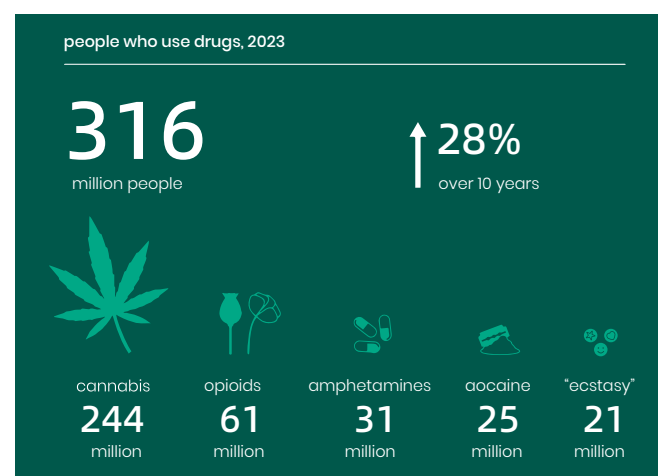

POCBio Drug of Abuse Screening Test

These tests offer rapid, reliable qualitative screening for one or more drugs and their metabolites in human urine.

Test Procedures


For Cassette

1



Remove the test cassette from the foil pouch.

2

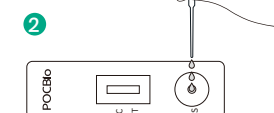

Dispense 3 drops of specimen vertically using the provided dropper to the "S" well.


3



Read the result between 5-10 minutes.


For Panel

1

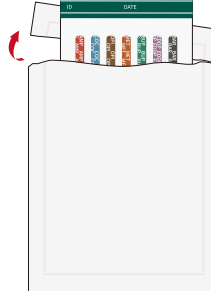

Open the foil pouch and remove the test device.

2

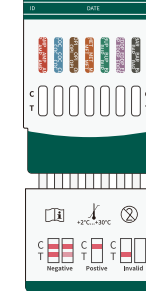


Remove the cap to expose the strips.

3



Collect the urine sample, then immerse the panel vertically into the sample for 10 seconds, ensuring the sample level is at the wavy line — not above or below.


4

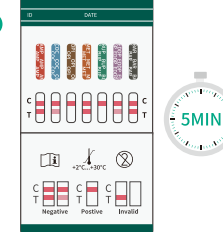

Place on a flat surface and read results at 5 minutes. Do not interpret after 10 minutes.


For Cup

1

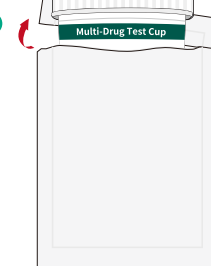


Remove the cup from the pouch.

2

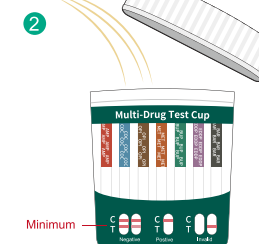


The donor urinates directly into the cup. Ensure the urine volume exceeds the minimum level.

3

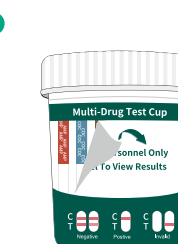

Securely re-cap the cup and place it on a flat surface. Peel off the label from left to right to view the results.

4

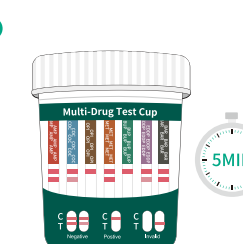

Read the results within 5 minutes. Do not interpret results after 10 minutes.


Result Interpretation



Negative



Positive

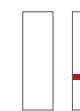

Invalid



Interpret adulteration strips between 3-5 minutes. Please refer to the color chart.